

TRU MARGIN™
RESTORATIVE ULTRASONIC INSTRUMENTS

Restorative Dentistry

TRU MARGIN™ ◆ DOUBLE DIAMOND ◆

Diamond Restorative Ultrasonic Tips


Use Tru Margin Restorative Tips for ultra precise prepping and fine finishing of the dentinal shoulder. The use of Tru Margin diamond coated ultrasonic instruments will enhance your preparations by creating a finer, gentler, more precise margin- a TRU margin.



To Order Toll Free:
(855) 737-9477 • (619) 470-2375
To Order Online: www.sdswiss.com
info@sdswiss.com • We Ship Worldwide

Patent Pending

SAN DIEGO
 **SWISS**

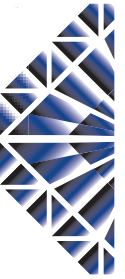

Double Diamond ◆ Water Port ◆ Contra Angle



FUNCTIONS AND USES

TRU MARGIN-1 (1mm x 1.4mm taper)

Radius preparation tip.

Use for preparation of the dentine (chamfer shape).

TRU MARGIN-2 (1.3mm x 1.7mm taper)

Radius preparation tip.

Use for preparation of the dentine and finishing the dentine (chamfer shape).

TRU MARGIN-3 (.6mm)

Flat shoulder preparation tip.

Use for preparation of the shoulder shaping.

TRU MARGIN-4 (1.1mm)

Flat shoulder preparation tip, larger at the distal end for shoulder shaping and finishing.

TRU MARGIN Ultrasonic Instruments are designed for a precise and accurate pre-prosthetics finish and a TRU Margin preparation for your crowns.

- **TRU MARGIN 1-2** Instruments - Radiused on the distal end of the instrument
- **TRU MARGIN 3-4** Instruments - Flat on the distal end of the instrument for accurate shoulder preparation